

Ordering Information

Part Number Description

Optiflex S Shoulder CPM

2029 OptiFlex S Shoulder CPM

Standard Accessories

34273 Control pendant

20661 Patient data chip card (1x)

Optional Accessories

20719 Grey Fleece Patient Kit (not shown)

86016 Retractable Casters

Technical Specifications

Mains Power: 100-240 V, 50/60 Hz, 33 VA

Weight: 25 kg

Dimensions*: 88 x 58 x 58 cm

Seat Height: 48 cm

Electrical Safety Class: Class 2, Type B

Safety Tests: IEC 60601-1, IEC 60601-1-2,
ANSI AAMI ES 60601-1,
CAN CSA 22.2 No. 60601-1-08

Range of Motion:

Adduction/Abduction: 0-30-175°

Internal/External Rotation: 90-0-90°

Flexion/Elevation: 0-30-175°

Horizontal Adduction/Abduction: 0-0-125°

Functions:

Motor Control: ON/OFF

Pauses: 0-30 seconds

Speed: 0-100%

Timer: 0-300 minutes

Synchronization: ON/OFF

Therapy Run Time: hours/minutes

* Dimensions are listed as L x W x H



190, Boul. St-Elzéar Ouest, Laval (Québec) H7L 3N3
T 450 669.8985 | 1 800 361.0877
F 450 669.9532 | 1 888 935.7001

www.mtm.ca

Optiflex® S Shoulder CPM



Not for sale in Europe

Features

OptiFlex S Continuous Passive Motion

Designed to provide anatomical motion for the shoulder, the OptiFlex S addresses a complete range of motion capabilities.


- Synchronized motors allow for customization of treatment parameters (optional to turn on/off)
- Patented universal left/right design for easy conversion
- Fully adjustable ergonomic chair that promotes shoulder stabilization
- Brushless DC motor
- Pendant control with chip card technology to store and assign patient information
- Quick and easy disassembly for transport and shipping
- 2 year warranty

Program Features

- Warm-up mode – the unit recognizes the start and stop angles and midpoint. From that point, motion is initiated 3 degrees from all angles until full ROM is achieved
- Reverse-on-load – unit will switch to the opposite direction when patient resistance exceeds the set value
- Patient data storage – assigns ROM, speed and treatment times per patient and notes patient set up positions
- Patient compliance meter – calculates the cumulative patient usage time
- Soft start and stop
- Multi-lingual pendant display



Adduction/
Abduction

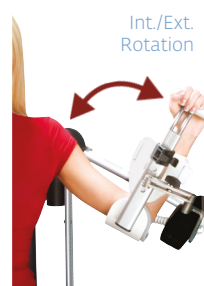

Int./Ext.
Rotation



Flexion/
Elevation



Physio
Movement